

INSTRUCTIONS FOR USE



■ HE0066NS Inactivated SARS-CoV-2 Whole Virus (Swab)

INTENDED USE

Inactivated SARS-CoV-2 Whole Virus (Swab) is intended for use as an external, non-viable control material to evaluate the performance of nucleic acid tests that detect SARS-CoV-2 virus. This product has no qualitative or quantitative assigned value. This control material is nonautomated and not intended to be used for screening, monitoring, or diagnosis. This control is not intended for any specific patient population or specimen.

SUMMARY AND PRINCIPLES

Inactivated SARS-CoV-2 Whole Virus (Swab) contains Severe Acute Respiratory Syndrome Coronavirus 2 isolate USA/WA1/2020. This control can be used to monitor the extraction, amplification and detection of SARS-CoV-2 molecular testing assays. Routine use of quality controls monitor test variation, lot to lot test kit performance, operator performance, and aid in identifying random or systemic error.

COMPOSITION

Inactivated SARS-CoV-2 Whole Virus (Swab) is comprised of cultured and inactivated Severe Acute Respiratory Syndrome Coronavirus 2 isolate USA/WA1/2020, and human A549 cells. The SARS-CoV-2 virus has been inactivated using chemical and heat treatments, and the A549 cells have been inactivated using thermal treatment. The viral particles are prepared in a buffered solution with materials of plant and animal origin, preservatives and stabilizers. The product consists of six individually packaged lyophilized swabs.

WARNINGS AND PRECAUTIONS

- For In Vitro Diagnostic use only.
- Do not open foil pouch until ready to use.
- This product must be handled in accordance with Biosafety Level 1 practices as described in the United States Department of Health and Human Services Centers for Disease Control and Prevention (CDC) and National Institutes of Health (NIH), Biosafety in Microbiological and Biomedical Laboratories, or other equivalent guidelines. Although this product has been inactivated; there is no known test or inactivation method that can assure that this product will not transmit infection.
- Wear proper personal protective equipment.
- Refer to the Safety Data Sheet (SDS) for more detailed information. The SDS can be located on the Microbiologics website at www.microbiologics.com or by contacting Technical Support at 320.229.7045 or U.S. Toll Free 1.866.286.6691.
- This product does not contain any hazardous substances listed in 67/548/EEC or listed in 1272/2008/EC.
- This product is not made with natural rubber latex.
- Customer is responsible for ensuring compatibility of the control with the assay or protocol in use.



MATERIALS REQUIRED BUT NOT PROVIDED

- Nucleic acid extraction kit and assay
- Instrumentation for detection
- Rehydration buffer such as nuclease-free water, sample preparation reagent, or transport medium as required by assay to be performed
- Pipettors capable of delivering 0.5-1000µl volumes
- Nuclease-free aerosol barrier pipette tips
- Vortex (if applicable)

INSTRUCTIONS FOR USE

1. Read assay package insert, instructions for use, or applicable lab protocol. Some instruments and assays are equipped with special QC settings. In these instances, it may be necessary to use the special setting when using QC sets and panels.
2. Tear open pouch at notch. Remove swab from pouch and place into a vial of hydrating fluid.
 - a. Minimum hydration volume of 300ul is required for complete hydration of the swab. Please note, the swab will absorb some of the hydration fluid.
 - b. Note: The swab may be hydrated in liquid or processed directly following diagnostic assay instructions for use.
3. Break the swab by snapping the shaft - lift the swab a few millimeters from the bottom of the vial and pushing against the rim to break it.
4. Recap the vial and vortex or shake vigorously for 10 seconds or until swab is hydrated.
 - a. Rehydrated material may be stored at 4°C and used up to 8 days after hydration. Mix before use.
5. Use the appropriate volume for the assay being performed and follow laboratory protocols or manufacturer instructions for processing a sample.
6. Note: Each swab is intended as a single use test. Dilutions may be performed and used immediately.

STORAGE AND EXPIRATION

Store the Inactivated SARS-CoV-2 Whole Virus (Swab) at 2°C - 25°C in the original packaging up to the indicated expiration date. After opening the foil pouch use the swab immediately.

Inactivated SARS-CoV-2 Whole Virus (Swab) should not be used if:

- Stored improperly
- There is evidence of excessive exposure to heat or moisture
- The expiration date has passed

LIMITATIONS

This product may not be suitable for use with all kits and procedures. Customer is responsible for verifying the performance of this product with their chosen instrumentation and assay(s).

MICROBIOLOGICAL STATE

This product was prepared using suitable inactivation methods. While the product has been tested for innocuity, universal laboratory precautions are recommended, and material should be treated as though it was a viable specimen.

KEY OF SYMBOLS

 LOT	Batch Code (Lot)		Manufacturer
 REF	Catalog Number		Temperature Limitation
	Caution, Consult Accompanying Documents		Use By
	Contains sufficient for <n> tests		Refer to Instructions for Use
 IVD	In Vitro Medical Device		Telephone Number
 EC REP	Authorized Representative in the European Community		CE Mark
		 QTY	Quantity

PRODUCT WARRANTY

- These products are warranted to meet the specifications and performance printed and illustrated in product inserts, instructions, and supportive literature.
- The warranty, expressed or implied, is limited when:
 - The procedures employed in the laboratory are contrary to printed and illustrated directions and instructions
 - The products are employed for applications other than the intended use cited in product inserts, instructions, and supportive literature

NOTICE TO PURCHASERS

The purchase of this product allows the purchaser to use it for Research and Quality Control. No general patents or other license of any kind other than this specific right of use from purchase is granted hereby. No other rights are conveyed expressly, by implication or by estoppel to any other patents. Furthermore, no rights for resale are conferred with the purchase of this product.

The Microbiologics logo and Helix Elite™ are registered trademarks of Microbiologics, Inc.

WEBSITE

Visit our website, www.microbiologics.com, for current technical information and product availability.

ACKNOWLEDGEMENTS



Microbiologics, Inc.

200 Cooper Avenue North
St. Cloud, MN 56303 USA

Customer Service

Tel. 320-253-1640
U.S. Toll Free 800-599-BUGS (2847)
Email. info@microbiologics.com

Technical Support

Tel. 320-229-7045
U.S. Toll Free 866-286-6691
Email. techsupport@microbiologics.com
www.microbiologics.com



Medimark® Europe

11, rue Emile Zola B.P. 2332
38033 Grenoble Cedex 2, France

Tel. 33 (0)4 76 86 43 22


Fax. 33 (0)4 76 17 19 82

Email. info@medimark-europe.com

ILLUSTRATED INSTRUCTIONS


Each kit consists of 6 individually packaged lyophilized swabs.

1

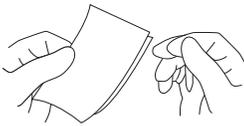

Read assay package insert, instructions for use, or applicable lab protocol. Some instruments and assays are equipped with special QC settings. In these instances, it may be necessary to use the special setting when using QC sets and panels.

2

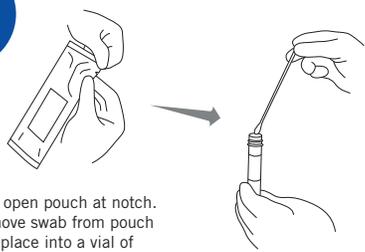

Tear open pouch at notch. Remove swab from pouch and place into a vial of hydrating fluid.

Minimum hydration volume of 300ul is required for complete hydration of the swab. Please note, the swab will absorb some of the hydration fluid.


Note: The swab may be hydrated in liquid or processed directly following diagnostic assay instructions for use.

3



Break the swab by snapping the shaft - lift the swab a few millimeters from the bottom of the vial and pushing against the rim to break it.

4

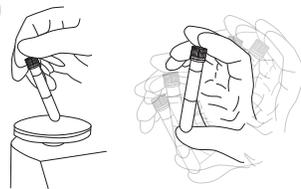


Recap the vial and vortex or shake vigorously for 10 seconds or until swab is hydrated.

Rehydrated material may be stored at 4°C and used up to 8 days after hydration. Mix before use.



5

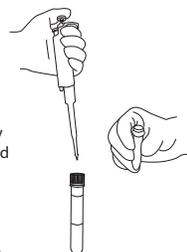

Use the appropriate volume for the assay being performed and follow laboratory protocols or manufacturer instructions for processing a sample.

6

Note: Each swab is intended as a single use test. Dilutions may be performed and used immediately.